


 Jiangsu Vedkang Medical Science & Technology Co., Ltd	Documentation no.	VDK/JS-CE-35-05
	Instructions for Use	English
	Effective Date: 2018.05.23	Rev. A/0

Polyp Trap

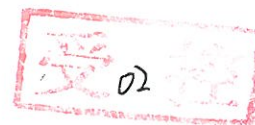

Instructions for Use of Polyp Trap

Originator/Date: 戴治平 2018.05.22

Reviewer/Date: 曹芳 2018.05.22

Approver/Date: 王华 2018.05.23

Effective Date: 2018.05.23

Status of document change control

SN.	Approved Date	Effective Date	Previous Version	Change No.	Change content
1	2018.05.23	2018.05.23	NA	20180523001	New product, so we use the A/0 version.


Instructions for Use of Polyp Trap

【Product Name】 Polyp Trap

【Structure】 Main parts: connecting tube, lid, strainer chamber, vessel and etc.

【Specification】 JH-PT-05

【Feature】 The product is made of PP, PVC, PC, etc. Integrative design of five chambers, independent number for each chamber, convenient for signing and marking, transparent shell, convenient for observing the placing situation of specimens.

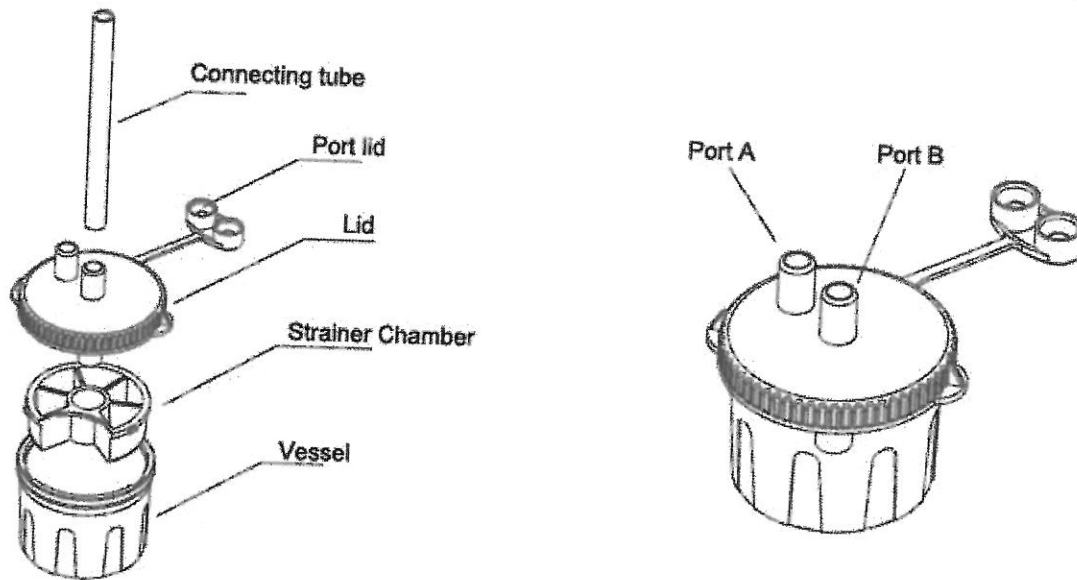

【Operations】

I

1. Take the polyp trap out from the package; check if any damage or the lid is closed tightly.
2. Connect Port A to the endoscope, and connect Port B to the suction system. Connecting tube may be used for transition if necessary.
3. Rotate to make the arrow on lid point to the number of the chamber, and then start suctioning. Record trapping place according to chamber number, and then fill related documents.
4. After specimen trapping, remove the connection tube between Port A and Port B.
5. Rotate to make the arrow on lid point to the Neutral position to pour out the waste fluid within the polyp trap.
6. Fill the polyp trap with formaldehyde solution and seal the Port A and Port B tightly with the port caps to avoid any volatilization and leakage.
7. Send out the polyp trap for inspection.

II

Dispose the product as per recognized hospital standards.

【Cautions】

1. Do not use this equipment when the package is damaged.
2. Do not use this product if the shelf-life is over.
3. Single use only, dispose after use and reuse is forbidden.
4. Prepare a spare Polyp Trap before surgery.
5. Connect the Polyp Trap, the suction source and the sampling channel firmly. Fail of doing this may affect the specimen collecting.
6. Fill the polyp trap with formaldehyde solution and seal the Port A and Port B tightly with the port caps to avoid any volatilization and leakage.
7. Please keep the products in a dry environment.
8. If the inner package of Polyp Trap has been damaged, do not use it.

【Package & Storage】

Inner packing uses Dupond bag; middle packing uses paper box; outer packing uses corrugated carton. The storage shall be in the room of well ventilation, the relative humidity is not higher than 80%, and has non-corrosive gas.

【Symbol】

- | | | |
|-----|---|---|
| 1. |  | Single use only |
| 2. |  | Non-sterilized product |
| 3. |  | Use by (expiration date) |
| 4. |  | Production lot number |
| 5. |  | Refer to instructions |
| 6. |  | CE mark |
| 7. |  | Manufacturing date |
| 8. |  | Do not use when package is damaged |
| 9. |  | Reference number |
| 10. |  | Manufacturer |
| 11. |  | Authorized representative in the European Union |

【Sterilization style】 EO sterilization

【Sterilization valid period】 Two years

【Production license number】 Jiangsu Food & Medicine Production No 2003-0066

【Name of EU representative】 Prolinx GmbH

【Address of EU representative】 Brehmstr. 56, 40239, Duesseldorf, Germany

【Telephone of EU representative】 0049 211 3105 4698

【Fax of EU representative】 0049 2131 4051968-9

【E-mail of EU representative】 med@eulinx.eu

【Name of manufacturer】 Jiangsu Vedkang Medical Science & Technology Co., Ltd

【Address of manufacturer】 No. 52, Guoxiang Road, Wujin Economic Development Zone, Changzhou, Jiangsu, China

【Telephone of manufacturer】 +86 519 69877751

【Fax of manufacturer】 +86 519 69877753

【Postcode of manufacturer】 213149

【E-mail of manufacturer】 sales@vedkang.com

【Manufacturing Date】: See Label

【Expiration Date】: See Label

【Version No】 A/0

【Revise Date】: May 22, 2018